

INTENDED PURPOSE

Commode wheelchairs are required when using or getting to the toilet is not possible due to the disability. Commode wheelchairs can be used as room chairs. Commode wheelchairs can be pushed over standard toilet bowls.



Marketed by Sanrai Med India!
Office Tower B, Elante Office No.506,
Plot No.178-178A, Industrial Area Phase 1,
Chandigarh: 160002

Customer Support: +91 7341103371
Email: customersupport@sanraimed.com
Website: www.sanraimed.com

WHEELED COMMUNE TRS 200 XXL



Wheeled Commode TRS 200—**your ultimate solution for mobility and convenience** in personal care. Designed with comfort in mind, this versatile commode **features smooth-rolling wheels, ensuring easy transport and use.** Ideal for those seeking independence and ease, the TRS 200 redefines **comfort and dignity in daily living.**

FEATURES

- ✓ Wide and stable wheeled commode for maximum requirements
- ✓ Wide and stable wheeled commode for maximum requirements
- ✓ Stable, folding, and removable footrests
- ✓ Fold-down armrests simplify transfer
- ✓ Parking brake on the rear casters
- ✓ Seat and cover can be folded up
- ✓ Continuously adjustable footrests
- ✓ With receptacle as standard
- ✓ Extremely simple assembly and cleaning



- | | |
|--------------------------------|----------------------------------|
| 1. Push handle | 7. Double-function parking brake |
| 2. Backrest | 8. Toilet bucket with cover |
| 3. Fastening clip for backrest | 9. Footrests |
| 4. Seat | 10. Heel strap |
| 5. Seat cover | 11. Footrest locking mechanism |
| 6. Armrest locking mechanism | 12. Armrest |

ACCESSORIES



Seat for TRS 200 XXL



Seat cover for TRS 200 XXL

TECHNICAL DATA


Product	TRS 200 XXL
Overall length	96 cm
Overall width	72 cm
Overall height	98 cm
Seat width	56 cm
Seat depth	42 cm
Seat height	52 cm
Height between floor and seat	44 cm
Weight	15.3 kg
Max. load	200 kg
Frame material	Welded steel tubing, powder-coated
Seat / armrest cover material	PP (polypropylene) / PP (polyurethane)

ASSEMBLY

- ✓ The wheeled commode can be easily assembled. Place the preassembled back unit on the base frame and secure it with the fastening clip (photo 2).


Locking and unlocking the armrests

- ✓ Pull both armrests upwards until they audibly engage (photo 3).
- ✓ They will automatically click into place. Ensure that they successfully click into place on both sides; in this case, the locking pin (1) must be visible (photo 4).
- ✓ To unlock the armrests, pull the locking knob back (2) and swing the armrest downwards (photo 4).

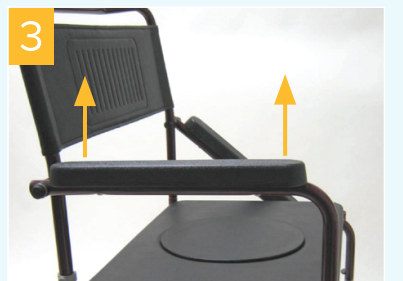


Double-function parking brake

- ✓ To apply the brake, press the lever downwards (photo 5). At the same time, the directional position of the wheel is blocked. To release the brake, press the head of the lever downwards (photo 6).

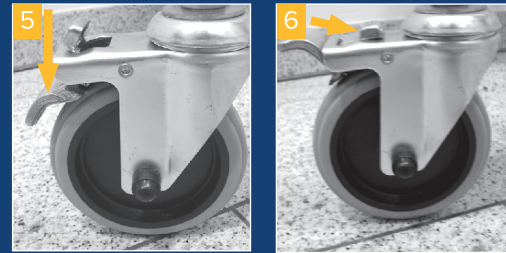

Mounting the footrests

- ✓ Hang the provided footrests with the two holes on the respective retaining pins as shown (photo 7+8)

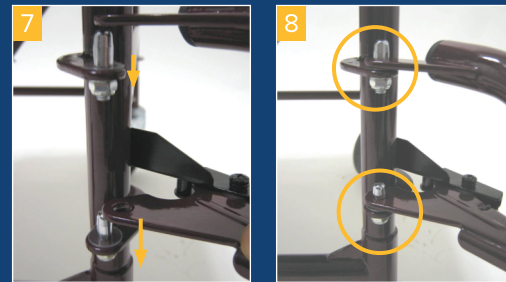

Locking and unlocking the footrests

- ✓ To lock the footrests, swivel them to the center until they audibly click into place. To unlock them, press the lever inwards as shown and move the footrests to the side (photo 9)



Toilet bucket

- ✓ Always hang and remove the toilet bucket from above (also with the cover on). The toilet bucket snaps into place and the toilet seat can be easily lifted off for cleaning (photo 10).



Height adjustment of the footrests

- ✓ To be able to adjust the height of the footrests, first loosen the clamping lever (photo 11+12). Then press both pins and adjust the footrests to the desired height (photo 13). Lock the clamping lever once again (photo 14).

IMPORTANT: Never use footrests to assist with standing up!



- ✓ The product should not be exposed to temperatures above 60°C for long periods of time (for example, use in a sauna, etc.)
- ✓ In the event of defects, contact your specialist dealer from whom you purchased this product.

- ✓ Do not use any high - dose chlorine and vinegar-based cleaners!
- Check the condition of the tires regularly. They are “fail-safe”, meaning that you do not have to check air pressure or fill them with air. Always check the tires for tears or other damage, however.

SAFETY INSTRUCTIONS

- ✓ Read these instructions for use carefully.
- ✓ All components should be checked for damage and a secure fit prior to use.
- ✓ IMPORTANT: If you need to move in reverse over thresholds with a person sitting in the wheeled commode, there is a risk of the armrest lock coming loose.
- ✓ When parking a wheeled commode, always use the parking brakes.
- ✓ Avoid wet conditions – risk of slipping!
- ✓ Prevent unauthorized use, for example, by children.
- ✓ The wheeled commode is significantly less stable than a regular wheelchair, therefore avoid leaning to the side – risk of tipping!
- ✓ Clothing or body parts can be pinched when sitting down, standing up, inserting the toilet bucket, as well as when raising the armrests and locking them into place.
- ✓ If you need to move over thresholds, do so slowly and carefully, otherwise the wheels can become damaged!
- ✓ During transfer into/out of the wheeled commode, the footrests must be swiveled out of the way and the locking brakes must be applied.
- ✓ Never use the footrests as assistance when getting up – risk of tipping!
- ✓ Take note of the permitted maximum weight of the user.
- ✓ This information can be found in the technical data in these instructions for use as well as on the identification plate.

CLEANING | MAINTENANCE

- ✓ Clean regularly with a moist cloth!
- ✓ Do not use abrasive cleaners!
- ✓ Do not use any aggressive cleaning agents!
- ✓ To disinfect, a spray disinfectant can be used such as Descosept AF or Bacillol AF.
- ✓ Check the screws regularly and tighten, if necessary.